



## Capsule Endoscopy System



# LESS IS MORE



We learn. We listen. We create.

We have been learning from our past successes and frustrations.

We have been listening to our valued customers' expectations fulfilled or yet to be met.

Now, we have created the smartest capsule endoscopy system, OMOM HD.

#### Reading cannot be easier.

Integrated with the evolving Artificial Intelligence and the most advanced optical technology, the latest OMOM HD capsule endoscopy system aims to assist healthcare professionals with unprecedented ease of interpreting capsule endoscopy studies.

Less Reading Time, More Evident Details.




## OMOM HD

Your Ultimate Choice of Small Bowel Capsule Endoscopy

## Reading made Smart

#### The new face - Vue Smart

With major upgrade on patient management center and reading screen, the new *Vue Smart* software offers healthcare providers even more user-friendly interface and stress-free reading experience.

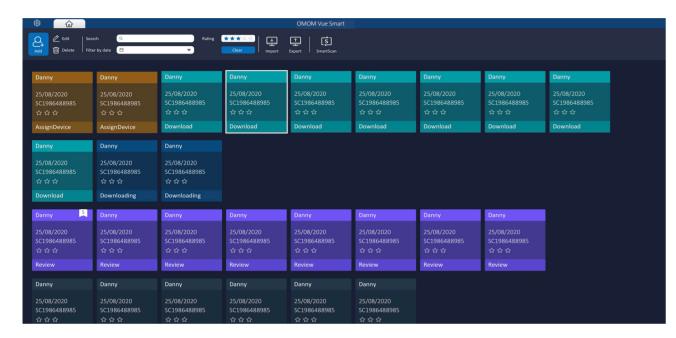




#### Al-Powered, Smarter Ever.

With our advanced AI algorithms, the *Vue Smart* is at the new technological frontier in redundancy deletion, lesion detection and classification that guarantees easier-than-ever interpretation as well as improved clinical outcomes.



Utilizing our proprietary algorithms, **SmartScan** processes the raw files by deleting up to 90% redundancy, identifying and tagging as many as 16 types of abnormalities. After *Smartscaning*, the *SmartView* and *SmartFinding* are ready for use.

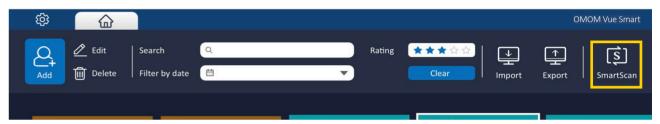




**SmartView** reading mode will only play the pictures picked up by *Smartscan*, assisting users in quickly navigating through the full video. Highlighted marks above the time bar show the images which can be displayed as a short video.





By capturing the most prominent characteristics of selected images with predictive results, **SmartFinding** helps doctors by displaying the abnormalities and suggesting descriptions. Images are displayed as page view.

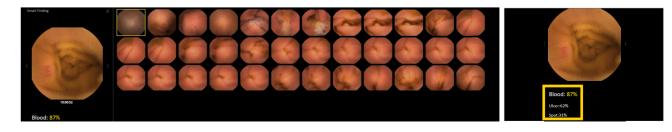

### Horizon Unbounded

Pushing the boundary on capsule hardware, our new OMOM HD delivers superior image quality powered by our new optical technology. With more evident details from higher resolution images, ultra-wide-angle view, SpeedSense technology and longer battery life, OMOM HD ushers in the era of simple but confident diagnosis.

#### **Industry-leading Features**



512\*512 Pixel

>300% improved resolution Unparalleled clarity



2-10 FPS

Enhanced sampling frequency. Upgrade by SpeedSense



172 Degrees

Ultra-wide viewing angle Improved tissue coverage



12 Hours

Increased capsule working hours Extended reach and duration



Compact capsule

Ergonomic design. Dimension 25.4mm\*11mm Weight 3g




Disposable belt cover

Disposable belt cover is now available as contact-free option.







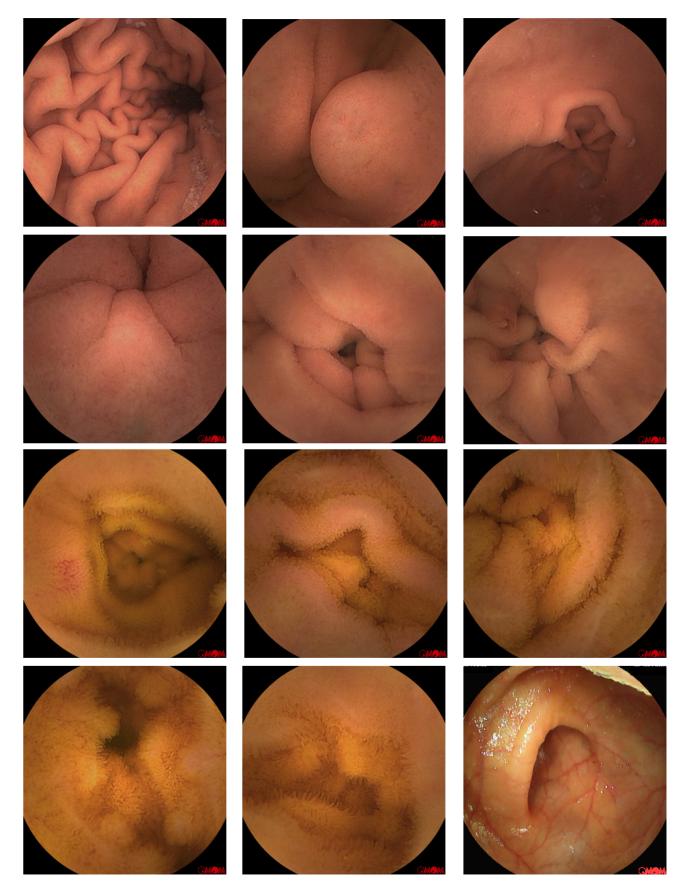





#### Chongqing Jinshan Science & Technology (Group) Co., Ltd.

Address: 18 Nishang Road, Yubei District, 401120 Chongqing, China

Tel: +86-23-86098151 Fax: +86-23-86098777 Email: international@jinshangroup.com

www.jinshangroup.com/en

CE Registration No: HD 60146955 0001 Specifications subject to changes without prior notice

